




Clinical Systems

Selectra Pro M

THE SMARTEST CHOICE
FOR YOUR GROWING
LABORATORY NEEDS

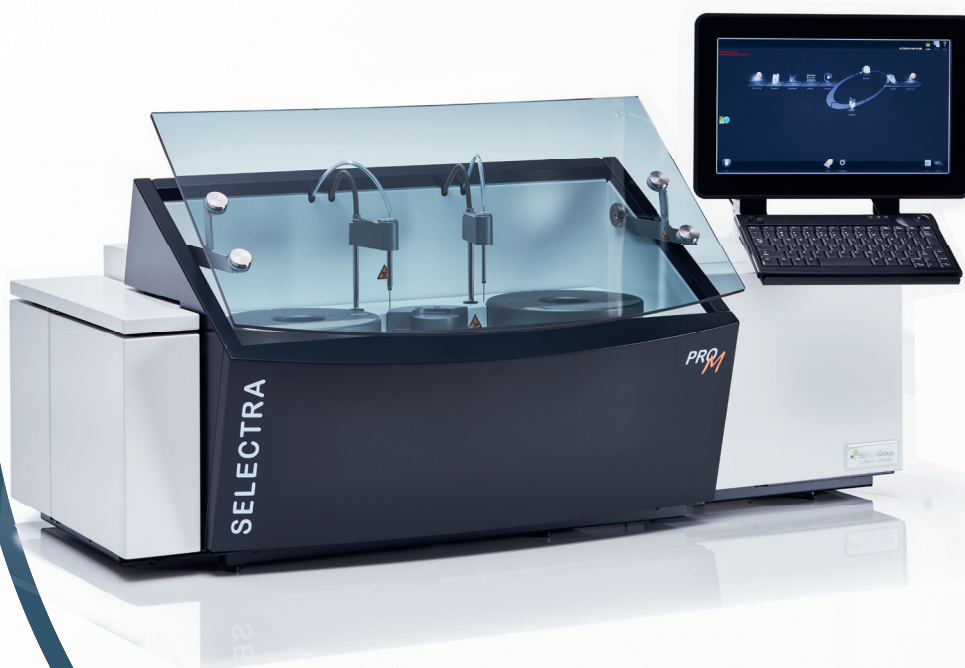

Selectra Solutions



SPECIFICATIONS

THROUGHPUT

- Typically 180 tests/hour
- Up to 266 ISE tests/hour (Dry ISE optional)

REAGENT HANDLING

- Refrigerated reagent rotor controlled with Peltier elements
- Cooled to 10°C+/-4°C at normal laboratory conditions
- 32 reagent positions for 10 mL, 25 mL and 50 mL reagent bottles
- Reagent ID and automated programming of assays, calibrators and controls via barcode
- All reagent positions can be assigned as R1, R2 and R3
- Pre-heated needle with level detection, collision protection and integrated mixer

SAMPLE HANDLING

- Sample rotor containing:
 - Outer segment with 50 barcode readable positions
 - Inner segment with 12 auxiliary positions
 - All positions fit 13x75 mm primary and secondary tubes and pediatric cups
 - All positions can be used for calibrators, controls, normal, pediatric and STAT samples
- Sample probe with level detection, integrated mixer and collision protection
- Programmable dilution ratios 1 : 5 up to 1 : 200 in one step increments with 3 possible diluents

PIPETTING SYSTEM

- 1000 µL reagent syringe:
 - R 1 volume 110 - 400 µL
 - R 2 volume 0 - 180 µL
 - R 3 volume 0 - 180 µL
 - Programmable in 1 µL steps
- 100 µL sample syringe:
 - Sample volume 1 - 30 µL
 - Programmable in 0.1 µL steps

CUVETTE ROTOR

- Cost effective, semi-disposable cuvette rotor with 48 cuvettes, path length 7 mm
- >10,000 tests per rotor
- Minimum measuring volume 220 µL
- Measuring temperature 37°C, controlled by Peltier elements

LIGHT SOURCE

- Quartz-iodine lamp 12V-20W

WAVELENGTH RANGE

- 340 - 800 nm
- Optical unit with 8 position filter wheel
- Automatic wavelength selection
- 340, 405, 505, 546, 578, 620, 660, 700 nm standard installed
- Other wavelengths available on request

PHOTOMETRIC RANGE

- -0.1 to 3.0 Absorbance
- Resolution 0.001 Abs

ANALYTICAL MODES (SINGLE, DUAL AND TRIPLE REAGENT SYSTEM)

- Kinetic measurement with linearity check
- Mono- and bichromatic end point measurement with or without bichromatic reagent blank and/or sample blank correction
- Two point measurement; with or without slope blank
- Graphic plot of all measuring points
- Predilution, post-dilution and automatic reflex dilution as needed
- Non-linear calibration curves
- Prozone check for immunology tests
- Cut-off declaration
- Calculated tests

QUALITY CONTROL

- Up to 15 different controls can be defined, 3 per test
- Westgard rules
- Levey-Jennings plots
- Quality control statistics

WATER CONSUMPTION

- ~500 mL per hour max, continuous operation

STANDARDS AND REGULATIONS

- CE - IVDD
- USA FDA 510(k)
- CB
- UL

DIMENSIONS & WEIGHT

- 125 cm (50 in) x 75 cm (30 in) x 62 cm (25 in) (W x H x D)
- 95 kg (210 lbs)

INTERFACE

- State of the art Host-Query interface available
- Host: RS232 or Ethernet (TCP/IP) through LIS-2A protocol
- Hand held CCD barcode reader used for reagent identification and automated programming of assays, controls, and calibrators

INSTALLATION CONDITIONS

- Temperature: 15 - 32 °C (59 - 90°F)
- Humidity: 15 - 85% RH
- Altitude: up to 2000 m
- Plumbing: no dedicated system water or drain required
- Electrical: Voltage: 100 - 240 Vac; Frequency: 50 / 60 Hz; Power (max): 400 VA

INTEGRATED PC

- Touch screen navigation
- Operating System: MS Windows™ Embedded

OPTIONS

ISE MODULE

- Patented Solid State Dry Electrode Technology
- Indirect measurement
- Dilution 1:14
- Measures Sodium, Potassium, Chloride and Bicarbonate

POSITIVE SAMPLE IDENTIFICATION

- Positive Sample Identification (PSID) via integrated barcode reader
- Reads all popular formats including Codes 39, 128, 11, 93, 4, CODABAR and Interleaved 2/5

PRINTER

- Printer supported by MS Windows™

PROACTIVE MAINTENANCE KIT

- Complete parts kit for annual preventive maintenance



WORLDWIDE OFFICES

Please contact your sales representative for terms, conditions and product availability in your country.

Headquarters

Australia
BELUX
Brazil
France
Italy

+33 1 41 45 07 10

1800 815 098
+32 92 820 531
+55 27 3025 1415
+33 4 83 36 10 82
+39 011 97 61 91

New Zealand
Serbia
The Netherlands
UK
United States

0800 555 611
+381 11 2467119
+31 313 430 581
+44 1442 869320
+1 800 453 2725



Order code: 6003-404



ELITechGroup

EMPOWERING IVD

www.elitechgroup.com
info@elitechgroup.com